

ScientekTM

Hospital and Laboratory Equipment

CAGE & BOTTLE WASHER



SW 5000 Series

The SW 5000 series washers from Scientek make for the perfect fit to your washing criteria. The flexibility and versatility of the interchangeable header system, allows you to process a full range of caging and glassware.

The simple yet well-proven design, coupled with 'state of the art' programming makes the Scientek SW 5000 series washers the right choice of washer for your facility.



**SW 5000 Control Panel
with PV 600 user interface**

- **The powerful and effective rotating spray washing arms, together with the versatile cycle programming, make this the winning washing combination.**
- **Simple touch control user interface**
- **Full display cycle status and diagnostics**
- **Pre-programmed 10 cycle selection**
- **Event, time and temperature displayed**
- **Modem connection for factory upgrades**




**Model SW 5200
c/w CVH-52
Combination Header**


Scientek's exceptional features include:

- All stainless steel type 304 welded construction on tank and pipework
- Sanitary gaskets and bolted mid seam for ease of installation and access
- Powerful all stainless steel close coupled pump
- Self cleaning incline debris strainer
- Guaranteed 180° F final rinse treatment
- Non-proprietary parts
- Safety door hardware and emergency controls
- Pneumatic control valves standard
- Drop door loading platform
- Proven effective overlap rotary spray washing
- Powerful impinging water pressure sprays ware clean
- Allen Bradley controls & program hardware

SW 5000 Series – Ware Washing Accessories:

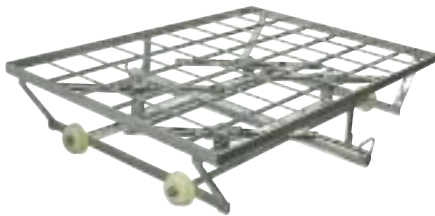

BRH-55
Bottom Rotary Header
(flat deck washing platform)

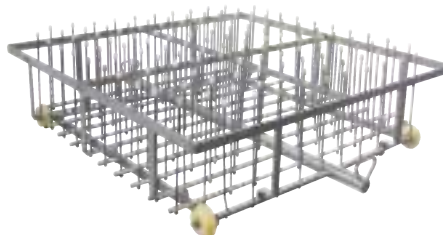

SH-55-140
Spindle Header
(140 position bottle washing)



TT-55
Transfer Table
(mobile header storage)



SMR-55
Vertical Header
(used between
cage racks)



SMR-55
Cage Rack
(two per load with
vertical header)



BB-36
36 Bottle Capacity Basket
(multiple basket options
available)

SW5000 CAGE & BOTTLE WASHER SERIES

Table of washer cabinet models for the size you need.

| MODEL # | INTERIOR CABINET SIZE | EXTERIOR CABINET SIZE | PUMP MOTOR |
|-----------|--|--|------------|
| SW 5200 | 47"W x 36"D x 40"H 1194 x 914 x 1016 | 69"W x 44 1/2"D x 80"H 1753 x 1130 x 2032 | 3.0 HP |
| SW 5200PT | 47"W x 36"D x 40"H 1194 x 914 x 1016 | 69"W x 46"D x 80"H 1753 x 1168 x 2032 | 3.0 HP |
| SW 5500 | 47"W x 54"D x 40"H 1194 x 1372 x 1016 | 69"W x 61"D x 80"H 1753 x 1549 x 2032 | 7.5 Hp |
| SW 5500PT | 47"W x 54"D x 40"H 1194 x 1372 x 1016 | 69"W x 62 1/2"D x 80"H 1753 x 1588 x 2032 | 7.5 Hp |

- **Note: Certain options can increase exterior sizes**
- **Scientek will provide mechanical service and installation drawings on request**
- **Please consult factory for further details on how we can help you.**

Optional washer features & accessories allow for personal customization:

- **Higher output 7 1/2HP pump (std. SW 5500)**
- **Power operated guillotine door/s**
- **Electric heating in place of standard steam heating**
- **Automatic cold water drain discharge cooling**
- **120° F temperature rise to process supply water**
- **Direct DI or pure water final rinse**
- **Recessed wall or barrier isolation**
- **Cycle data printer**
- **Data system reporting**
- **Direct quick connect accessory washing (std. SW 5500)**
- **CGMP upgrades to washer**
- **Chemical dispensing and monitoring system**
- **Devapourmatic system – self contained (no exterior exhaust)**
- **Power operated exhaust vent**
- **Dryer system – self contained (no exterior exhaust)**
- **SMR-55 - Cage Washing System (2 racks & vertical header)**
- **BRH-55 - Bottom Rotary Spray Header (flat deck washing platform)**
- **CVH-55 - Combination Header (bottles and cages)**
- **SH-55 - Spindle Headers (multiple combinations and sizes)**
- **TT-55 - Transfer Header Accessory Table**

For additional product information please contact the factory or the Scientek representative in your area.



**11151 Bridgeport Road
Richmond, BC, Canada V6X 1T3
Tel: (604) 273-9094
Fax: (604) 273-1262
E-mail: sales@scientek.net
www.scientek.net**